

## Troubleshooting guide

| Problem                          | Causes   | Comments and suggestions   |
|----------------------------------|--|--|
|                                  | Upside down during transportation may cause beads to stick with sealing film                           | Spin down the 'Cartridge' by hand before open it.  |
| Physical damage of the kit       | Insufficient reagent volume due to evaporation or lifting of the sealing film                          | If reagent volumes are insufficient, extraction efficiency may decrease. Do not use it and contact customer service immediately.   |
|                                  | Broken 'Cartridge'   | Do not use it and contact customer service immediately.  |
| Inappropriate specimen condition | Specimen condition is not favorable due to inappropriate storage condition (ex. coagulation)           | Perform vortex and pipetting of specimen. If sample is still coagulated, add a bit of PBS or distilled water and vortex again.   |
|                                  | Specimen condition is abnormal due to inappropriate storage condition (ex. stored in room temperature) | Increase sample volume up to 300 µl, if extraction efficiency is low.  |
| Abnormal extraction              | Too much beads left in Elution buffer  | If the total nucleic acid density is within the normal range, proceed with the eluted solution. In the case of low total nucleic acid density, transfer the eluted solution to a 1.5 ml tube and centrifuge before use.  |
|                                  | Eluted total nucleic acid should not appear transparent or sticky                                      | Refer to 'Inappropriate specimen condition' part of trouble shooting if specimen condition is abnormal, perform extraction again. If the specimen condition and total nucleic acid density are normal, proceed with the extracted total nucleic acid. In cases where the specimen condition is normal but the total nucleic acid density is too sticky, transfer the eluted solution to 1.5 ml tube and dilute it with distilled water before use. |

## Warnings and precautions

- Intended for in vitro diagnostics.
- Intended for professional use only.
- Read and follow the Instruction for Use (IFU) before using the product.
- Use extracted nucleic acid as soon as possible, if long-term storage is needed, store it below -70 °C.
- Be cautious of contaminants such as microorganisms after opening the product.
- Be sure to wear personal protective equipment such as gloves and goggles when using this product and wash hands after handling specimens and reagents.
- Avoid contamination with DNase or RNase during product use.
- Store the product at the specified storage temperature and do not use it beyond its expiration date.
- Read and follow the IFU for the nucleic acid extraction device (AllEx® Mini Automated Nucleic Acid Extraction System) used with this product.
- Do not dispose of reagents from this product with bleach or acidic substances, as they contain irritants.
- This product is a single use and should not be reused.

\* A notice to the user that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or the patient is established.

## Storage conditions

- Temperature : 15 °C to 25 °C
- Relative humidity : 20 % to 80 %



GeneAll Bldg., 303-7 Dongnam-ro, Songpa-gu, Seoul, 05729, Republic of Korea  
E-mail : sales@geneall.com  
Tel : 82-2-407-0096 Fax : 82-2-407-0779



A-1201~A-1204, Hanam Techno Valley U1 Center, 947,  
Hanam-daero, Hanam-si, Gyeonggi-do, 12982, Republic of Korea

## GENEALL BIOTECHNOLOGY CO., LTD.

©2025 GeneAll Biotechnology, All rights reserved

| Symbol | Used for                     | Symbol | Used for            |
|--------|------------------------------|--------|---------------------|
|        | Batch number                 |        | Manufacturer        |
|        | Catalogue number             |        | Do not re-use       |
|        | Consult instructions for use |        | Date of manufacture |
|        | Caution                      |        | Expiry date         |
|        | Temperature limit            |        |                     |

2025.03

www.geneall.com

Ver. 2.0

GeneAll®

Store between 15 °C and 25 °C  
Expiration date: Refer to the product box label

Allex® Mini  
Fecal DNA/RNA Kit

## Description

The AllEx® Mini Fecal DNA/RNA Kit is a specialized nucleic acid extraction reagent designed for use with the AllEx® Mini Automated Nucleic Acid Extraction System (AllEx® Mini). The Kit protocols are optimized to extract DNA and RNA using AllEx® Mini Cartridge, which can extract up to 12 sample at once. Each cartridge comes preloaded with essential reagents, enabling dependable and efficient nucleic acid extraction while minimizing user effort.

## Intended purpose

AllEx® Mini Fecal DNA/RNA Kit provides fast and easy methods for the purification of DNA and RNA from various fecal samples using AllEx® Mini Automated Nucleic Acid Extraction System.

The extracted nucleic acid is of excellent quality and suitable for most downstream applications including PCR, qPCR, RT-PCR, or any other downstream application.

## Kit contents

| Components                      | 976-048      |
|---------------------------------|--------------|
| Number of preparation           | 48 preps/kit |
| Reagents pre-filled cartridge   | 6 pks        |
| Buffer FL 70 ml                 | 1 ea         |
| 2 ml Glass bead tube (50 ea/pk) | 1 pk         |
| Proteinase K 33 mg *            | 1 ea         |
| PK-Storage buffer 2 ml          | 1 ea         |



Cartridge




Cassette

\* Prepare Proteinase K (20 mg/ml) solution according to the instruction before use.

\* After reconstitution of Proteinase K store at 2 °C to 8 °C.

\* For long-term storage, it is recommended to store it below -20 °C.

## Brief workflow

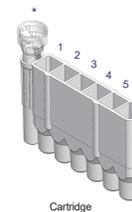

Cartridge

- Lysis buffer
- Washing buffer I
- Washing buffer II
- Washing buffer III + bead

[Note] If particles are not visible in well 4, shake down the cartridge to dislodge particles that may have adhered to the seal material before removing the seal.

- Elution buffer
- Washing buffer IV
- Strip

- AllEx® Mini Cartridge is designed for low throughput, allowing flexible processing of up to 12 individual samples.
- The Strip pocket has a Strip that keeps the magnet safe from magnetic beads while making sure samples and reagents mix well.
- Six consecutive wells are aligned horizontally, each containing specific reagents for extraction.
- The first well contains lysis buffer that breaks down cell membranes and releases DNA and RNA.
- To prevent damage, the magnetic beads are placed in the fourth well until the extraction process begins, at which point the AllEx® Mini transfers them to the first well. The second, third, fourth, and sixth well contain washing buffer I, II, III, and IV, which remove unwanted cellular components and residual buffers. Finally, the elution buffer in the fifth well separates DNA and RNA from the magnetic beads, completing the extraction process.

## ■ Protocol


The kit includes two protocols, allowing selection based on the specific extraction purpose.

| Protocol                     | Uses and purpose  |
|------------------------------|---|
| P1 Protocol<br>(10 min 51 s) | <ul style="list-style-type: none"> <li>• Rapid, efficient and PCR-compatible nucleic acid extraction</li> </ul>     |
| P2 Protocol<br>(20 min 47 s) | <ul style="list-style-type: none"> <li>• High-quality nucleic acid extraction for NGS-grade applications</li> </ul> |

### A. Fecal


1. Transfer 0.2 g of fecal sample to a 2 ml Glass bead tube.
2. Add 1 ml of Buffer FL and 25 µl of Proteinase K Solution (20 mg/ml). Pulse vortex the mixture at least 1 min.
3. Incubate at 65 °C for 5 min.
4. Centrifuge at 13,000 rpm (≥ 10,000 xg) for 5 min at room temperature.
5. Transfer 200 µl of the supernatant to 1<sup>st</sup> well of Kit.

## ■ Preparation of Cartridge

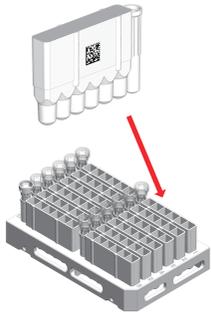

1. Insert up to 12 Cartridges in the desired quantity, starting with position number 1 indicated on the Cassette.
3. Remove the sealing film.
4. Follow the protocol for starting sample type.


## ■ System run

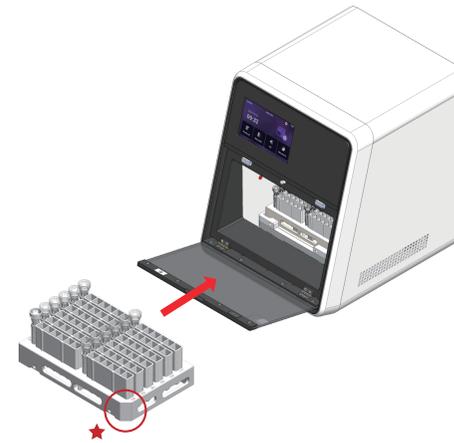

1. Load the Cassette with the inserted Cartridges into the AllEx® Mini, ensuring the diagonally cut edge faces the system door.




2. Tap the play icon to scan the barcode on the Cartridge.

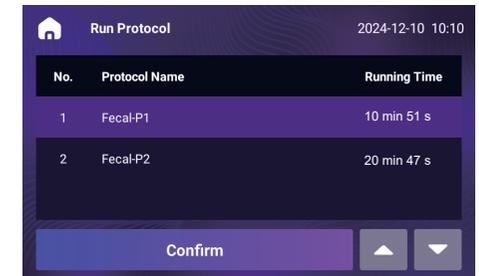

3. The built-in Barcode reader reads the Cartridge barcode, selects protocol and displays the protocol lists. Select the protocol from the list based on uses and purposes.
4. Tap 'Confirm' to start extraction.

## ■ End of run

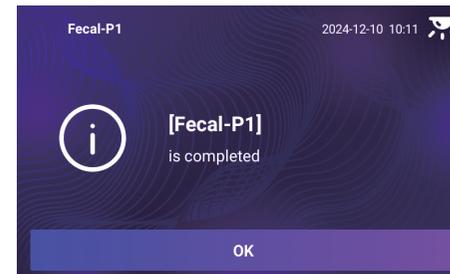

1. Once extraction is complete, tap 'OK' to return to the main menu.
2. Open door and check whether the Strip is correctly placed in the Strip pocket.
3. Remove the Cassette from AllEx® Mini.
4. Regular UV sterilization eliminates bacteria and virus, minimizing internal pollution within the system. Before/after the extraction, close the door and proceed with UV sterilization as a decontamination activity. The sterilization will continue for 10 min.